

A decorative vertical strip on the left side of the page, featuring a complex pattern of white, overlapping, curved lines that create a sense of depth and movement against a light blue background.

RELIANCE ADHESIVES

LIGHT BOND™

LIGHT CURE ORTHODONTIC ADHESIVE



CONTENTS OF KIT:

- 20 G PASTE
- (2) 3 CC SEALANT RESIN
- 23 G LIQUID ETCHING AGENT
- MIXING PAD
- SPATULAS
- BRUSHES (SPONGE PELLETS MAY BE SUBSTITUTED FOR BRUSHES)

PRODUCT FEATURES:

- **AVAILABLE WITH OR WITHOUT FLUORIDE**
LIGHT BOND™ PASTE AND SEALANT ARE FORMULATED WITH A PATENTED MONOMER THAT CONTAINS FLUORIDE. IT IS AVAILABLE WITHOUT FLUORIDE IN PUSH SYRINGE ONLY.
- **UNSURPASSED BOND STRENGTH**
THE LIGHT BOND™ PASTE IS A HIGHLY FILLED MATERIAL THAT DELIVERS MAXIMUM BOND STRENGTH IN ALL SITUATIONS WITH ALL TYPES OF APPLIANCES
- **SMOOTH TACKY PASTE VISCOSITY**
THE VISCOSITY OF LIGHT BOND™ PASTE IS VERY TACKY WHICH PREVENTS BRACKET FLOTATION. IF A THINNER PASTE VISCOSITY IS DESIRED FOR BONDING LINGUAL RETAINERS AND ACRYLIC APPLIANCES, LIGHT BOND™ THIN OR MEDIUM PASTE IS AVAILABLE. SPECIFY WHEN ORDERING.
- **TOTAL PASTE POLYMERIZATION UNDER ALL BRACKET BASES**
LIGHT BOND™ IS COMPRISED OF A CATALYST AND FILLER THAT REFLECTS THE LIGHT TO INSURE COMPLETE POLYMERIZATION UNDER ANY APPLIANCE: METAL, CERAMIC, SAPPHIRE OR PLASTIC. ONCE THE MATERIAL HAS BEEN PROPERLY CURED, AN ARCH WIRE CAN BE PLACED IMMEDIATELY.
- **EASE OF REMOVAL**
LIGHT BOND™ IS A SMALL PARTICLE GLASS FILLED MATERIAL DESIGNED SPECIFICALLY FOR ORTHODONTIC USE THAT ALLOWS MAXIMUM PENETRATION OF FILLER INTO THE FINEST SCREEN MESH BASE AND REMOVAL WITHOUT ENAMEL SCARRING.
- **SUGGESTED USES**
BRACKET ADHESIVE

LIGHT BOND KITS

LIGHT BOND™ FILLED SEALANT MAY BE SUBSTITUTED FOR STANDARD LIGHT BOND™ SEALANT. NON-FLUORIDE RELEASE FORMULA AVAILABLE UPON REQUEST.

LIGHT BOND™ PASTE IS AVAILABLE IN REGULAR, MEDIUM AND THIN VISCOSITIES.

LIGHT BOND™ CAN BE USED WITH ANY CURED SEALANT AND IS AVAILABLE IN KIT FORM WITH ASSURE® UNIVERSAL BONDING RESIN (HYDROPHILIC) OR LIGHT BOND™ SEALANT.

LIGHT BOND™ IS AVAILABLE IN KIT FORM WITH PRO SEAL® OR L.E.D. PRO SEAL® ENAMEL PROTECTIVE SEALANT FOR AN ADDITIONAL CHARGE.

BLUE SHADE AVAILABLE UPON REQUEST.

LIGHT BOND™ KIT WITH PASTE IN PUSH SYRINGES WITH FLUORIDE	AD-1113K
LIGHT BOND™ ECONOMY KIT WITH PASTE IN LONG NOSE TIPS WITH FLUORIDE	AD-1115K
LIGHT BOND™ ECONOMY KIT WITH PASTE IN PUSH SYRINGES NO FLUORIDE	AD-1118K
LIGHT BOND™ KIT WITH PASTE IN PUSH SYRINGE WITH FLUORIDE AND ASSURE SEALANT	AD-1114K
LIGHT BOND™ ECONOMY KIT WITH PASTE IN PRELOADED TIPS WITH FLUORIDE	AD-1126K
LIGHT BOND™ KIT WITH PASTE IN SCREW SYRINGES WITH FLUORIDE	AD-1116K
LIGHT BOND™ BLUE KIT WITH PASTE IN PUSH SYRINGES WITH FLUORIDE	AD-1105K
LIGHT BOND™ BLUE ECONOMY KIT WITH PASTE IN PRELOADED TIPS WITH FLUORIDE	AD-1107K
LIGHT BOND™ KIT MEDIUM PASTE IN PRELOADED TIPS WITH FLUORIDE	AD-1112K
LIGHT BOND™ MEDIUM PASTE KIT WITH PASTE IN PUSH SYRINGES WITH FLUORIDE	AD-1111K
LIGHT BOND™ THIN PASTE KIT WITH PASTE IN PUSH SYRINGES WITH FLUORIDE	AD-1127K
LIGHT BOND™ THIN PASTE ECONOMY KIT WITH PASTE IN PRELOADED TIPS WITH FLUORIDE	AD-1130K

INDIVIDUAL LIGHT BOND™ COMPONENTS

5 G LIGHT BOND™ PASTE IN PUSH SYRINGE, WITH FLUORIDE	AD-1120
5 G LIGHT BOND™ PASTE IN PUSH SYRINGE, WITHOUT FLUORIDE	AD-1124
5 G LIGHT BOND™ PASTE IN PRELOADED TIPS, WITH FLUORIDE	AD-1123
5 G LIGHT BOND™ PASTE IN LONG NOSE TIPS, WITH FLUORIDE	AD-1119
5 G LIGHT BOND™ PASTE IN SCREW SYRINGE, WITH FLUORIDE	AD-1122
5 G LIGHT BOND™ MEDIUM PASTE IN PRELOADED TIPS, WITH FLUORIDE	AD-1117
5 G LIGHT BOND™ MEDIUM PASTE IN PUSH SYRINGE, WITH FLUORIDE	AD-1121
5 G LIGHT BOND™ THIN PASTE IN PUSH SYRINGE, WITH FLUORIDE	AD-1128
5 G LIGHT BOND™ THIN PASTE IN PRELOADED TIPS, WITH FLUORIDE	AD-1129
3 CC LIGHT BOND™ SEALANT WITH FLUORIDE	AD-1103
3 CC LIGHT BOND™ SEALANT WITHOUT FLUORIDE	AD-1125
7 CC LIGHT BOND™ SEALANT WITH FLUORIDE	AD-1108
3 CC LIGHT BOND™ FILLED SEALANT WITH FLUORIDE	AD-1109
7 CC LIGHT BOND™ FILLED SEALANT WITH FLUORIDE	AD-1110

INDIVIDUAL LIGHT BOND™ BLUE COMPONENTS

5 G LIGHT BOND™ BLUE PASTE IN PUSH SYRINGE, WITH FLUORIDE	AD-1104
5 G LIGHT BOND™ BLUE PASTE IN TIPS, WITH FLUORIDE	AD-1106

PAD LOCK®

LIGHT CURE FLUORESCING ADHESIVE



CONTENTS OF KIT:

- 18 G PASTE
- (2) 3 CC SEALANT RESIN
- 9 G LIQUID ETCHING AGENT (LIQUID OR GEL)
- MIXING PAD
- SPATULAS
- BRUSHES (SPONGE PELLETS MAY BE SUBSTITUTED FOR BRUSHES)

REFRIGERATION IS RECOMMENDED

PRODUCT FEATURES:

- **BONDS CHEMICALLY TO METAL**
PAD LOCK® WAS DEVELOPED TO INCREASE ADHESION TO EVEN THE FINEST SCREEN MESH BASES. THE PASTE CONTAINS A MONOMER THAT BONDS CHEMICALLY AS WELL AS MECHANICALLY TO METAL BRACKET BASES IMPROVING ADHESION.
- **FLUORESCES**
PAD LOCK® PASTE WILL FLUORESCENE WHEN ILLUMINATED WITH ANY ULTRAVIOLET BLACK LIGHT. THIS FEATURE ALLOWS THE OPERATION TO DISCLOSE FLASH NOT REMOVED AT THE TIME OF BRACKET PLACEMENT AND ADHESIVE REMNANTS LEFT AT DEBONDING.
- **SINGLE PASTE, LIGHT ACTIVATED**
- **COMPATIBLE WITH ALL SEALANTS AND SELF-ETCHING PRIMERS**
- **AVAILABLE WITH OR WITHOUT FLUORIDE**
- **SUGGESTED USES**
BRACKET ADHESIVE.

PAD LOCK® KITS

ASSURE® UNIVERSAL BONDING RESIN OR PRO SEAL® ENAMEL PROTECTIVE SEALANT MAY BE SUBSTITUTED FOR STANDARD LIGHT BOND™ SEALANT FOR AN ADDITIONAL CHARGE.

PAD LOCK® KIT IN PRELOADED TIPS WITH FLUORIDE AND ASSURE SEALANT **AD-1132K**

PAD LOCK® KIT IN PRELOADED TIPS WITHOUT FLUORIDE, AND ASSURE SEALANT **AD-1134K**

PAD LOCK® KIT IN PRELOADED TIPS WITHOUT FLUORIDE **AD-1133K**

PAD LOCK® KIT IN PRELOADED TIPS WITH FLUORIDE **AD-1135K**

PAD LOCK® KIT IN PUSH SYRINGES WITHOUT FLUORIDE AND ASSURE SEALANT **AD-1137K**

PAD LOCK® KIT IN PUSH SYRINGES WITH FLUORIDE AND ASSURE SEALANT **AD-1138K**

PAD LOCK® KIT IN PUSH SYRINGE WITH FLUORIDE **AD-1139K**

PAD LOCK® KIT IN PUSH SYRINGE WITHOUT FLUORIDE **AD-1136K**

INDIVIDUAL PAD LOCK® COMPONENTS

4.5 G PAD LOCK® PASTE IN PUSH SYRINGE WITH FLUORIDE **AD-1200**

4.5 G PAD LOCK® PASTE IN PUSH SYRINGE WITHOUT FLUORIDE **AD-1140**

4.5 G PAD LOCK® PASTE IN TIPS WITH FLUORIDE **AD-1141**

4.5 G PAD LOCK® PASTE IN TIPS WITHOUT FLUORIDE **AD-1201**

3 CC LIGHT BOND™ SEALANT WITH FLUORIDE **AD-1103**

ALSO AVAILABLE

6 ML L.E.D. PRO SEAL® LIGHT CURE FLUORESCING FILLED SEALANT **AD-1148**

6 ML PRO SEAL® LIGHT CURE FLUORESCING FILLED SEALANT **AD-1149**

ASSURE™

HYDROPHILIC LIGHT CURE ORTHODONTIC ADHESIVE

CONTAINS FLUORIDE



CONTENTS OF KIT:

- 20 G PASTE IN SCREW SYRINGES
- (1) 6 CC ASSURE® UNIVERSAL BONDING RESIN
- 9 G ETCHANT (LIQUID OR GEL)
- MIXING PADS
- SPATULAS
- BRUSHES (SPONGE PELLETS MAY BE SUBSTITUTED FOR BRUSHES)

REFRIGERATION REQUIRED

PRODUCT FEATURES:

- **HYDROPHILIC**
ASSURE® PASTE IS A LIGHT CURE PASTE THAT WILL BOND TO SLIGHTLY CONTAMINATED ENAMEL. ASSURE® PASTE IS THE IDEAL VISCOSITY TO PREVENT BRACKET FLOTATION AND CONTAINS FLUORIDE.
- **BONDS TO ANY ENAMEL TYPE**
ASSURE® UNIVERSAL BONDING RESIN WILL BOND TO ANY ENAMEL TYPE, WET OR DRY, NORMAL OR ATYPICAL, INCLUDING DECIDUOUS, BLEACHED AND DENTIN. IN ADDITION, ASSURE® WILL BOND TO COMPOSITE AND METAL SURFACES WITHOUT ANY ADDITIONAL PRIMER.
- **NO LIGHT CURING REQUIRED**
CAN BE USED UNDER ANY CHEMICAL, LIGHT OR DUAL CURE ADHESIVE, REGARDLESS OF MANUFACTURER.
- **SUGGESTED USES**
BRACKET ADHESIVE, WET FIELD BONDING.

ASSURE® KIT

COMPLETE ASSURE® KIT

AD-1090K

INDIVIDUAL ASSURE® COMPONENTS

6 CC ASSURE® UNIVERSAL BONDING RESIN

AD-1092

5 G ASSURE® PASTE IN SCREW SYRINGE

AD-1091

2 ML ASSURE® UNIVERSAL BONDING RESIN

AD-1093

GoTo™

ADHESIVE PASTE



GoTo SYRINGE

ADHESIVE PASTE IN 4G PUSH SYRINGE

AD-1225



GoTo TIPS

ADHESIVE PASTE IN 10 PRELOADED TIPS

AD-1226

ASSURE PLUS™

LIGHT CURE ORTHODONTIC ADHESIVE



ASSURE® PLUS

6 CC ASSURE® PLUS UNIVERSAL BONDING RESIN

AD-1215

ASSURE™ Z

PRODUCT FEATURES:

ASSURE UNIVERSAL BONDING RESIN WILL INCREASE ADHESION TO NORMAL, ATYPICAL, DRY OR SLIGHTLY CONTAMINATED ENAMEL. THIS UNIVERSAL BONDING AGENT BONDS TO METAL RESTORATIONS WITHOUT ADDITIONAL PRIMERS. NO DRYING OR LIGHT CURING REQUIRED AND CAN BE USED WITH ANY LIGHT OR CHEMICAL CURE SYSTEM.

ASSURE UNIVERSAL BONDING RESIN WILL BOND TO THE FOLLOWING SURFACES WITHOUT ADDITIONAL PRIMERS:

- FLUOROSED
- DECIDUOUS
- BLEACHED OR WET ENAMEL
- DENTIN
- METAL SURFACES
- COMPOSITE RESTORATIONS

ASSURE UNIVERSAL BONDING RESIN IS ENHANCE ADHESION BOOSTER AND BONDING RESIN IN ONE STEP!

PRODUCT FEATURES:

NEW ASSURE PLUS ALL SURFACE BONDING RESIN WILL INCREASE ADHESION TO NORMAL, ATYPICAL, DRY OR SLIGHTLY CONTAMINATED ENAMEL AND NOW ZIRCONIA AND ACRYLIC TEMPS WITHOUT ADDITIONAL PRIMERS. THIS UNIVERSAL BONDING AGENT BONDS TO METAL RESTORATIONS AND NOW PORCELAIN WITHOUT HYDROFLUORIC ACID. CAN BE USED WITH ANY LIGHT OR CHEMICAL CURE SYSTEM.

ASSURE PLUS ALL SURFACE BONDING RESIN WILL BOND TO THE FOLLOWING SURFACES:

- FLUOROSED
- DECIDUOUS
- BLEACHED OR WET NORMAL AND ATYPICAL ENAMEL
- DENTIN / CEMENTUM
- GOLD, AMALGAM, STAINLESS STEEL
- COMPOSITE RESTORATIONS
- ACRYLIC TEMPORARY / PONTIC TEETH
- PORCELAIN



ASSURE®

6 CC ASSURE® UNIVERSAL BONDING RESIN

AD-1091

ENHANCE™

ADHESION BOOSTER FOR ATYPICAL OR CONTAMINATED ENAMEL



PRODUCT FEATURES:

ENHANCE™ ADHESION BOOSTER WILL INCREASE THE ADHESION OF ANY CHEMICAL, LIGHT OR DUAL CURE SYSTEM TO THE FOLLOWING SURFACES:

FLUOROSED, HYPO-CALCIFIED OR DECIDUOUS ENAMEL, ACRYLIC, BLEACHED OR WET ENAMEL, DENTIN, METAL SURFACES, COMPOSITE RESTORATIONS AND CONTAMINATED ENAMEL.

INDIVIDUAL ENHANCE® COMPONENTS

3 CC ENHANCE™ ADHESION BOOSTER PART A **AD-1012**

3 CC ENHANCE™ ADHESION BOOSTER PART B **AD-1013**

6 CC ENHANCE™ ADHESION BOOSTER PART A **AD-1014**

6 CC ENHANCE™ ADHESION BOOSTER PART B **AD-1015**

ENHANCE™ L.C.

SINGLE COMPONENT ENHANCE LIGHT CURE



PRODUCT FEATURES:

- PROVIDES ALL THE BENEFITS OF CHEMICAL CURE ENHANCE™ WITHOUT THE MIXING. MUST BE USED WITH LIGHT CURE SEALANT. THIS SINGLE LIQUID (NO MIXING) IS APPLIED TO THE ETCHED ENAMEL SURFACE THEN LIGHTLY AIR DRIED. LIGHT BOND™ SEALANT IS PLACED DIRECTLY ONTO THE ENHANCE™ L.C. AND LIGHT CURED FOR 10 SECONDS (20 SECONDS REQUIRED IF PRO SEAL® IS USED).
- MOST CONVENIENT FOR LIGHT BOND™ AND PAD LOCK® USERS!

INDIVIDUAL ENHANCE® L.C. COMPONENT

6 CC ENHANCE™ LIGHT CURE **AD-1016**

PROTECTIVE EYEWEAR

MORE AND MORE ORTHODONTISTS ARE USING LIGHT CURED ADHESIVES AND NEWER, MORE POWERFUL CURING LIGHTS. IT MAKES SENSE FOR YOU, YOUR ASSISTANTS AND PATIENTS TO PROTECT YOUR EYES FROM THE POWERFUL LIGHT SOURCES NOW BEING USED EVERYDAY IN YOUR CLINIC. RELIANCE HAS SOURCED A VERY HIGH QUALITY, YET VALUE PRICED PROTECTIVE EYEWEAR FROM UVEX. THEY ABSORB 99.9% OF HARMFUL UVA AND UVB RADIATION, AND ARE AVAILABLE IN TWO DESIGNS.



UVEX SKYPER® FOR USE WITHOUT PRESCRIPTION GLASSES

AD-1184



UVEX ULTRA SPEC® 2000 SIZED TO BE USED OVER YOUR PRESCRIPTION GLASSES

AD-1185

L.E.D. PRO SEAL® OR PRO SEAL® LIGHT CURE FLUORESCING FILLED SEALANT

- PRO SEAL® AND L.E.D. PRO SEAL® ARE THE FIRST SEALANTS THAT WILL COMPLETELY SET WITHOUT AN OXYGEN INHIBITED LAYER. THIS CREATES A SMOOTH HARD SURFACE THAT PROTECTS THE ENAMEL.
- PRO SEAL® AND L.E.D. PRO SEAL® ARE HIGHLY FILLED, NO MIX LIGHT CURE SEALANTS THAT RESIST TOOTHBRUSH ABRASION.
- PRO SEAL® AND L.E.D. PRO SEAL® ALSO CONTAIN A UV FLUORESCING AGENT FOR EASY MONITORING OF SEALANT COVERAGE.
- PRO SEAL® AND L.E.D. PRO SEAL® CAN BE USED UNDER ANY LIGHT CURE, CHEMICAL CURE OR DUAL CURE PASTE SYSTEM. U.C.S.F. AND M.C.G. TESTS PROVE THE LONG TERM BENEFITS OF THIS ENAMEL PROTECTIVE SEALANT.



6 ML L.E.D. PRO SEAL®

AD-1148

6 ML PRO SEAL®

AD-1149

L.C.R.™

LIGHT CURE RETAINER ADHESIVE



CONTENTS OF KIT:

- 20 G PASTE
- (2) 3 CC SEALANT RESIN
- 23 G LIQUID ETCHING AGENT
- MIXING PADS
- SPATULAS
- BRUSHES (SPONGE PELLETS MAY BE SUBSTITUTED FOR BRUSHES)

PRODUCT FEATURES:

• LONGEVITY AND DURABILITY

L.C.R.™ LIGHT CURED BONDING PASTE IS A HYBRID COMPOSITE THAT CURES TO A SMOOTH GLASS FINISH FOR PATIENT COMFORT AND REDUCED WEAR. THIS HIGHLY FILLED LIGHT CURED RESIN PASTE IS AN EXCELLENT CHOICE WHEN LONGEVITY AND DURABILITY ARE REQUIRED. SIMPLY APPLY ACCORDING TO INSTRUCTIONS AND LIGHT CURE.

L.C.R.® KITS

L.C.R.™ KIT WITH PASTE IN PRELOADED TIPS

AD-1099K

L.C.R.™ KIT WITH PASTE IN PUSH SYRINGES

AD-1098K

INDIVIDUAL L.C.R.® COMPONENTS

5 G L.C.R.™ PASTE IN PUSH SYRINGE

AD-1101

5 G L.C.R.™ PASTE IN PRELOADED TIPS

AD-1102

1.5 G L.C.R.™ PASTE IN PUSH SYRINGE

AD-1100

FLOW TAIN™ AND FLOW TAIN™ L.V.

LIGHT CURE FLOWABLE MICROFIL COMPOSITE



SUGGESTED USES:

- LINGUAL RETAINER ADHESIVE
- ATTACHMENTS ON THERMOPLASTIC TOOTH ALIGNING PRODUCTS
- SPLINTING MATERIAL (E.G. FIBERS)
- INDIRECT BONDING ADHESIVE
- ARCH WIRE STOPS

PRODUCT FEATURES:

- **FLOW TAIN™ AND FLOW TAIN™ L.V.**
ARE FLOWABLE LIGHT CURED COMPOSITES. THESE MATERIALS ARE VERY SMOOTH WHEN CURED WHICH WILL REDUCE WEAR AND ENHANCE PATIENT COMFORT.
- **CURED COMPOSITE FLEXIBILITY**
FLOW TAIN™ AND FLOW TAIN™ L.V., WHEN CURED, HAVE A LOW MODULUS OF ELASTICITY AND ARE NOT AS BRITTLE AS CONVENTIONAL COMPOSITES.
- **FLOWABLE VISCOSITY**
EVEN THOUGH FLOW TAIN™ AND FLOW TAIN™ L.V. ARE FLOWABLE MATERIALS, THEIR THIXOTROPIC PROPERTIES PREVENT THEM FROM SLUMPING FOR BETTER CONTROL.
- **NEEDLE TIP DISPENSER**
USING THE TIPS PROVIDED, THESE MATERIALS ARE EASILY INJECTED ONTO LINGUAL SURFACES AND ACROSS FIBER SPLINTS.

INDIVIDUAL FLOW TAIN® COMPONENTS

1.5 G FLOW TAIN® WITH TIPS	AD-1095
1.5 G FLOW TAIN® LOW VISCOSITY® WITH TIPS	AD-1094
FLOW TAIN® .019 GAUGE TIPS	AD-1096
FLOW TAIN® LOW VISCOSITY .022 GAUGE TIPS	AD-1097

PERFECT-A-SMILE™

SINGLE LIGHT CURE PASTE



PRODUCT DESCRIPTION:

Perfect-a-Smile is a single light cure paste that comes in 5 shades and can be painted on any aligner material to mask a space or missing tooth. Perfect-a-Smile bonds to any clear thermoplastic aligner, no conditioner or preparation is needed!

INDIVIDUAL SHADES	
B-1	AD-1212
A-1	AD-1231
A-2	AD-1213
A-35	AD-1214
LIGHT	AD-1407
MEDIUM	AD-1409



B-1



A-2



A-3

ORTHODONTIC RETENTION

POST-TREATMENT RETAINING PRODUCTS

ORTHO FLEXTech - GOLD RETENTION SYSTEM



Ortho FlexTech is intended to be used for retaining tooth position after orthodontic treatment.

- Stabilizes orthodontic corrections
- Prevents rotational relapses
- Naturally forms to lingual arch curvature
- Less doctor time / less chair time
- Low failure rate due to "flex linkage"
- Improved patient comfort - lays flat against the teeth
- Lasting aesthetics due to 14 karat white gold

ORTH FLEXTech - DENTAL RETENTION

AD-1208

BOND-A-BRAID LINGUAL RETENTION WIRE



Permanent retention is greatly enhanced using a flattened, dead soft 8-braided wire. Intra-arch splinting with this thin ribbon arch prevents torque control problems that can occur using a round braided wire for these procedures. When bonded on the lingual of the upper incisors, a very thin rectangular wire can be adapted with little lower incisor interference.

The dead soft wire is easy to adapt and prevents inadvertent tooth movement associated with active force wires.

Bond-A-Braid is excellent for semi-permanent orthodontic splinting of the upper or lower incisors. It can be used to hold diastemas closed, maintain difficult extraction sites and solving various other tooth-splinting dilemmas. Excellent for splinting periodontally compromised teeth.

BOND-A-BRAID LINGUAL RETENTION SYSTEM

AD-1189

STAINLESS STEEL ORTHO-FLEXTech



130% tensile strength and 300% torsional strength compared to gold ortho-flextech!!!

STAINLESS STEEL ORTHO-FLEXTech

AD-1205

RETAINIUM WIRES



Single strand NICKEL FREE titanium flat metal ribbon with increased interproximal strength is ideal for semi permanent splinting and maintaining difficult extraction sites.

60" total: 10 x 6" strands

RETAINIUM LINGUAL RETENTION WIRES

AD-1210

BAND LOK™

TWO PASTE COMPOMER BAND CEMENT

CONTAINS FLUORIDE



CONTENTS OF KIT:

- 24 G PASTE A
- 24 G PASTE B
- MIXING PAD
- SPATULAS

PRODUCT FEATURES:

- **NO MESSY POWDER-LIQUIDS**
CONSISTENT MIX VISCOSITY EVERY TIME.
- **BONDS TO METAL**
BAND-LOK® WILL BOND CHEMICALLY TO STAINLESS STEEL BANDS AND CROWNS. REMOVAL IS ALL IN ONE PIECE WITH CEMENT IN THE APPLIANCE.
- **EXTENDED WORKING TIME**
FOUR MINUTE WORKING TIME AT ROOM TEMPERATURE WITHOUT FROZEN SLAB.
- **DUAL CURABLE**
BAND-LOK® IS A DUAL CURE CEMENT. IT CAN BE USED WITH OR WITHOUT DENTAL CURING LIGHT. AFTER FLASH IS CLEANED A 20-SECOND LIGHT ACTIVATION ALLOWS FULL PROTECTION FROM MOISTURE. WITHOUT LIGHT ACTIVATION A COMPLETE SET IS ACHIEVED IN 10 MINUTES.

AVAILABLE IN NATURAL OR BLUE SHADE

BAND LOK® KITS

BAND-LOK® KIT IN NATURAL SHADE

AD-1078K

BAND-LOK® KIT IN HIGH VISIBILITY BLUE SHADE

AD-1077K

INDIVIDUAL BAND LOK® COMPONENTS

6 G BAND-LOK® PASTE A IN SCREW SYRINGE

AD-1074

6 G BAND-LOK® PASTE B IN SCREW SYRINGE
(USED FOR BOTH NATURAL AND BLUE SHADES)

AD-1072

6 G BAND-LOK® PASTE A
BLUE SHADE IN SCREW SYRINGE

AD-1073

6 G BAND-LOK® PASTE DUAL CURE A&B SELF-
MIXING SYRINGE IN HIGH VISIBILITY BLUE SHADE

AD-1071

6 G BAND-LOK® PASTE DUAL CURE A&B
SELF-MIXING SYRINGE IN NATURAL SHADE

AD-1070

SELF MIXING DUAL SYRINGES FOR BAND-LOK® AND CROWN-LOK®

AVAILABLE WITH BAND-LOK®, BAND-LOK® BLUE, AND CROWN LOK® BLUE PASTE



CONTENTS OF KIT:

- 4 DUAL SYRINGES OF PASTE
- 60 MIXING TIPS

THESE SYRINGES MIX AND DISPENSE THE PASTE IN ONE EASY STEP. THE A + B SYRINGES ARE JOINED AND USE SPECIAL MIXING CANULA TIPS, WHICH ELIMINATES OPERATOR MIXING AND PROVIDE A CONSISTENT, AUTOMATIC MIX EVERY TIME.

SELF MIXING DUAL SYRINGE KITS

BAND-LOK® KIT CHEMICAL CURE A+B SELF-MIXING DUAL SYRINGE IN NATURAL SHADE

AD-1076K

BAND-LOK® KIT CHEMICAL CURE A+B SELF-MIXING DUAL SYRINGE IN HIGH VISIBILITY BLUE SHADE

AD-1075K

CROWN-LOK® KIT BLUE SHADE

AD-1035K

CROWN LOK™

LUTING CEMENT FOR BANDS AND METAL CROWNS



CONTENTS OF KIT:

- 24 G CATALYST PASTE A
- 24 G BASE PASTE B
- MIXING PADS
- 9 G ETCHING AGENT (LIQUID OR GEL)
- SPATULAS

REFRIGERATION REQUIRED

CROWN LOK™ IS A SELF CURING, TWO MINUTE WORKING TIME FLUORIDE RELEASING CEMENT. THIS VERSATILE CEMENT WILL BOND CHEMICALLY TO STAINLESS STEEL CROWNS AND BANDS, REDUCING THE CHANCE OF FAILURE AND FACILITATION CLEANUP AT THE CONCLUSION OF TREATMENT.

AVAILABLE IN NATURAL OR BLUE SHADE

CROWN LOK® KITS

CROWN-LOK® KIT IN NATURAL SHADE

AD-1038K

CROWN-LOK® KIT IN HIGH VISIBILITY BLUE SHADE

AD-1037K

INDIVIDUAL CROWN LOK® COMPONENTS

6 G CROWN-LOK® PASTE A + B IN SELF-MIXING DUAL SYRINGE

AD-1039

6 G CROWN-LOK® PASTE A IN SCREW SYRINGE

AD-1042

6 G CROWN-LOK® PASTE A BLUE SHADE IN SCREW SYRINGE

AD-1041

6 G CROWN-LOK® PASTE B IN SCREW SYRINGE

AD-1036

ULTRA BAND LOK™

SINGLE PASTE COMPOMER BAND CEMENT

ALSO AVAILABLE IN LUER LOC SYRINGE / **CONTAINS FLUORIDE**

PRODUCT FEATURES:

- **SINGLE PASTE SYSTEM**
NO WASTEFUL MIXING
- **LIGHT ACTIVATED**
NO WORKING TIME CONSTRAINTS
- **NO MIXING**
- **ALL PURPOSE CEMENT**
ULTRA BAND-LOK® CAN ALSO BE USED TO BOND LARGE ACRYLIC APPLIANCES AND AS AN OCCLUSAL BUILD-UP MATERIAL. WHEN BONDING THE APPLIANCES YOU MUST FIRST ETCH THE ENAMEL SURFACE.
- **BONDS TO METAL**
ULTRA BAND-LOK® WILL BOND CHEMICALLY TO STAINLESS STEEL BANDS. REMOVAL IS ALL IN ONE PIECE WITH CEMENT IN THE APPLIANCE

AVAILABLE IN NATURAL OR BLUE SHADE



CONTENTS OF LUER LOC ULTRA BAND-LOK KIT:

- 30 G PASTE IN LEUR LOC PUSH SYRINGES
- (60) 15 GAUGE LEUR LOC DISPENSING TIPS



ULTRA BAND LOK® KITS

ULTRA BAND-LOK® LUER LOC KIT WITH DISPENSING TIPS

AD-1151K

ULTRA BAND-LOK® LUER LOC KIT IN HIGH VISIBILITY BLUE SHADE WITH DISPENSING TIPS

AD-1150K

INDIVIDUAL ULTRA BAND LOK® COMPONENTS

5 G ULTRA BAND-LOK® SYRINGE

AD-1155

5 G ULTRA BAND-LOK® HIGH VISIBILITY BLUE SYRINGE

AD-1152

5 G ULTRA BAND-LOK® LUER LOC SYRINGE WITH DISPENSING TIPS

AD-1154

5 G ULTRA BAND-LOK® LUER LOC HIGH VISIBILITY BLUE SYRINGE WITH DISPENSING TIPS

AD-1153

(60) 15-GUAGE LUER LOC DISPENSING TIPS

AD-1131

GLASSLOK™

DUAL CORE GLASS IONOMER BAND CEMENT

CONTAINS FLUORIDE



CONTENTS OF STARTER KIT:

- 1 GLASS LOK LIQUID 15 G BOTTLE
- 1 GLASS LOK POWDER 15 G
- 1 CLEAR SCOOP
- 1 MIXING ACCESSORY PACK

CONTENTS OF ECONOMY KIT:

- 3 GLASS LOK LIQUID 25 G BOTTLE
- 1 GLASS LOK POWDER 100 G
- 1 CLEAR SCOOP
- 1 MIXING ACCESSORY PACK

PRODUCT FEATURES:

- SECURE BAND ADHESION
- HIGH COMPRESSIVE AND TENSILE STRENGTH
- BONDS IN WET OR DRY ENVIRONMENT
- RESISTS WASHOUT FROM UNDER BANDS
- CONTAINS FLOURIDE - PROTECTS ENAMEL
- BONDS CHEMICALLY TO METAL, REDUCING REMOVAL CLEANUP TIME
- DUAL CURE - CHEMICAL OR LIGHT ACTIVATED



GLASS LOK® KITS

GLASS LOK® STARTER KIT - NATURAL SHADE

AD-1384K

GLASS LOK® STARTER KIT - BLUE SHADE

AD-1414K

GLASS LOK® ECONOMY KIT - NATURAL SHADE

AD-1415K

GLASS LOK® ECONOMY KIT - BLUE SHADE

AD-1416K

INDIVIDUAL PRECEDENT® COMPONENTS

15 G GLASS LOK® LIQUID - BLUE SHADE

AD-1418

15 G GLASS LOK® LIQUID - NATURAL SHADE

AD-1417

25 G GLASS LOK® LIQUID - BLUE SHADE

AD-1420

25 G GLASS LOK® LIQUID - NATURAL SHADE

AD-1419

100 G GLASS LOK® POWDER

AD-1421

15 G GLASS LOK® POWDER

AD-1422

GEL ETCHANT IN SYRINGE



6.5 G GEL ETCH IN SYRINGES AND 10 TIPS

AD-1008

ETCHANT TIPS PKG OF 10

AD-1009

23 G ECONOMY ETCHANT LIQUID



23 G ECONOMY ETCHANT LIQUID

AD-1011

23 G ECONOMY ETCHANT LIQUID CHERRY

AD-1010

9 G ETCHANT LIQUID



9 G ETCHANT LIQUID

AD-1017

18 G ECONOMY ETCHANT GEL



18 G ECONOMY ETCHANT GEL

AD-1007

RELIANCE S.E.P.™

SELF ETCHING PRIMER



4 ML (MIXED ETCHANT AND PRIMER), 1 DISPENSER, 2 REFILLS, 3 MIXING WELLS

AD-1022K

SELF ETCHING PRIMER DISPENSER

AD-1021

SELF ETCHING PRIMER REFILL CARTRIDGE

AD-1025

RELIANCE S.E.P.™ SELF ETCHING PRIMER PROMO KIT INCLUDES: 1 DISPENSER, 1 REFILL, 1 MIXING WELL

AD-1024K

9 G ETCHANT GEL



9 G ETCHANT GEL

AD-1006

ETCHMASTER SMALL ABRASIVE TIPS - 25 pack

AD-1206

ETCHMASTER CO2 REFILL

AD-1207

8 CC PORCELAIN CONDITIONER

FOR BONDING TO CERAMIC CROWNS



8 CC PORCELAIN CONDITIONER

AD-1005

10 CC PLASTIC CONDITIONER

FOR BONDING TO COMPOSITE RESTORATIONS
AND ACRYLIC APPLIANCES



10 CC PLASTIC CONDITIONER

AD-1004

7 G PORC-ETCH™

PORCELAIN ETCHING AGENT WITH BARRIER GEL



7 G PORC-ETCH™ AND 3 CC BARRIER GEL

AD-1019

3 CC BARRIER GEL

AD-1020

8 ML METAL PRIMER

FOR BONDING TO GOLD, STAINLESS STEEL AND AMALGAM



8 ML METAL CONDITIONER

AD-1018

RESTORE™

DIAMOND POLISHING PASTE



2 OZ JAR

AD-1158

1ST + FINAL™

PRE-PACKAGED
ORTHODONTIC PUMICE



PRODUCT FEATURES:

- **PRE-PACKAGED**
200 STERILE CUPS FOR INDIVIDUAL PATIENTS USE. PRE-MIXED – NO MESSY MIXING WITH WATER. JUST OPEN AND APPLY WITH PROPHY CUP.
- **IMPROVED BOND STRENGTH**
UNIVERSITY TEST SHOWS THAT 1ST + FINAL® WILL PUMICE MORE EFFECTIVELY & RINSE CLEANER FOR INCREASED BOND STRENGTH.
- **NO ADDED CHEMICALS**
NO FLAVORING, FLUORIDES OR OILS WHICH COULD LEAVE A FILM THAT MAY INHIBIT OR WEAKEN THE BOND.
- **ENAMEL POLISHER**
1ST + FINAL® CAN BE USED AS THE FINAL POLISHING STEP WITH THE RELIANCE RENEW™ FINISHING SYSTEM.

1ST + FINAL®

1ST + FINAL® 200 STERILE CUPS

AD-1001

PREP STEP™

DISPOSABLE PROPHY ANGLE



PRODUCT FEATURES:

144 INDIVIDUALLY WRAPPED AND DISPOSABLE PREP STEP™ PROPHY ANGLES ARE MANUFACTURED TO RUN SMOOTH AND VIBRATION FREE. LATEX FREE, 4 WEB AND SOFT CUP ARE DESIGNED TO FIT A SLOW SPEED HAND PIECE. JUST THROW AWAY AFTER USE!

PREP STEP®

PREP STEP™ 144 PROPHY ANGLES

AD-1003

RENEW™ FINISHING SYSTEM BURS

THESE BURS ARE THE **FIRST STEP** IN THE RENEW™ FINISHING SYSTEM

RENEW™ BUR 118S

THE 118S IS A SPECIAL 18 FLUTED BUR DESIGN FOR CONTROLLED REMOVAL OF BULK BONDING ADHESIVE AFTER THE BRACKETS HAVE BEEN REMOVED. THESE NON-TAPERED FRICTION GRIP BURS CUT QUICKLY & DO NOT PROMOTE FISSURES.



RENEW™ BUR 118S

AD-1160

RENEW™ BUR 218

THE 218 BUR IS AN 18 FLUTED DESIGN THAT IS TAPERED TO REACH INTO THOSE TIGHT AREAS ALONG THE GINGIVAL & OCCLUSAL WHEN REMOVING RESIDUAL ADHESIVE. THE WIDE FLUTES SELF CLEAN.



RENEW™ BUR 218

AD-1161

RENEW™ BUR 118L

THE 118L IS A LONGER SHANK VERSION OF THE 118S WITH THE SAME EXCELLENT ADHESIVE REMOVAL PROPERTIES. THESE NON-TAPERED FRICTION GRIP BURS CUT QUICKLY & DO NOT PROMOTE FISSURES.



RENEW™ BUR 118L

AD-1159

RENEW™ LINGUAL BUR 815

THE 815 LINGUAL COMPOSITE BUR WITH ITS FRICTION GRIP, CUSTOM LENGTH AND CUTTING HEAD, MAKE THIS THE IDEAL CARBIDE BUR FOR LINGUAL COMPOSITE REMOVAL.



RENEW™ BUR 815

AD-1164

RENEW™ FINISHING SYSTEM POINTS #383

THE #383 POINT IS THE SECOND STEP IN THE RENEW™ FINISHING SYSTEM

THE RENEW™ FINISHING SYSTEM POINTS WILL REMOVE SMALL RESIDUAL AMOUNTS FROM THE TOOTH SURFACE AND WILL POLISH THE ENAMEL WITHOUT DAMAGE. IT'S RUBBERIZED, TAPERED SHAPE FACILITATES REMOVAL OF ADHESIVE CLOSE TO THE GINGIVAL. THE #383 POINT IS THE SECOND STEP IN THE RENEW™ FINISHING SYSTEM. THE #383 POINTS ARE AVAILABLE IN STRAIGHT HAND PIECE, LATCH OR FRICTION GRIP.

EACH PACK CONTAINS 12 POINTS.

STRAIGHT HAND PIECE



STRAIGHT HAND PIECE PACK OF 12

AD-1188

FRICTION GRIP



FRICTION GRIP PACK OF 12

AD-1186

LATCH



LATCH PACK OF 12

AD-1187

DIAMOND BUR 129


A MEDIUM ABRASION DIAMOND BUR RECOMMENDED FOR ROUGHENING AND PREPARING PLASTIC, COMPOSITE RESTORATIONS AND METAL SURFACES. RELIANCE NOW OFFERS AN EXCELLENT QUALITY MEDIUM DIAMOND BUR AT A DOWN TO EARTH PRICE.



DIAMOND BUR 129, SOLD INDIVIDUALLY

AD-1163

THE COMPLETE ORTHODONTIC BONDING REFERENCE GUIDE

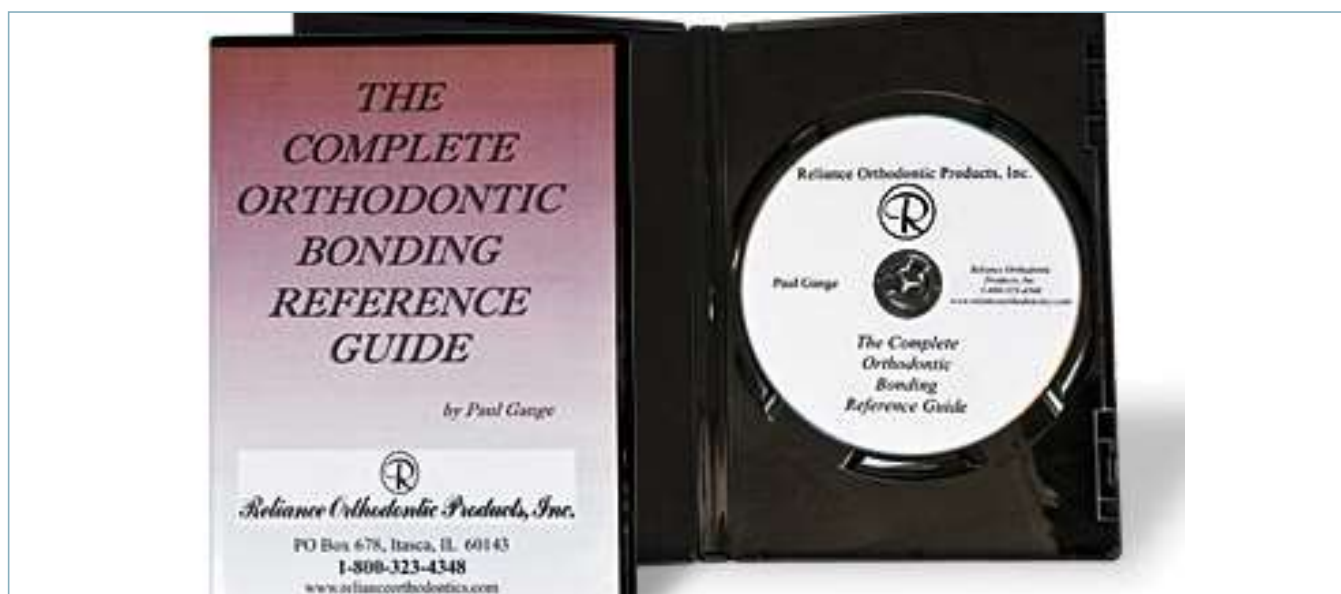

REFERENCE GUIDE INDEX:

- SECTION 1 **TOOTH PREPARATION**
- SECTION 2 **DIRECT BONDING ADHESIVES**
 - NO MIX
 - TWO PASTE
 - LIGHT CURE
 - DUAL CURE
- SECTION 3 **BONDING TROUBLESHOOTING**
- SECTION 4 **BONDING TO COMPOSITE RESTORATIONS,
METAL SURFACES AND VENEERS**
- SECTION 5 **BONDING IN A WET FIELD**
- SECTION 6 **BONDING LINGUAL RETAINERS**
 - USE OF FLOWABLE MICROFIL
 - OCCUSAL BUILDUP
- SECTION 7 **BAND CEMENTATION**
- SECTION 8 **CROWN CEMENTATION**
- SECTION 9 **BONDING LARGE ACRYLIC APPLIANCES**
- SECTION 10 **DEBONDING & ADHESIVE REMOVAL**
- SECTION 11 **INDIRECT BONDING**

PRODUCT FEATURES:

DEVELOPED BY PAUL GANGE, RESPECTED EXPERT ON ORTHODONTIC BONDING, THE COMPLETE ORTHODONTIC BONDING REFERENCE GUIDE COVERS EVERY BONDING TOPIC AND INCLUDES TOOTH PREPARATION, BONDING, DEBONDING, LINGUAL RETENTION, INDIRECT BONDING AND EVERYTHING IN BETWEEN! IDEAL FOR TRAINING CHAIR SIDE STAFF.

EXCELLENT REFERENCE FOR DOCTOR AND EXPERIENCED ASSISTANTS. COMPREHENSIVE, SUCCINCT AND EASY TO USE, THIS STATE-OF-THE-ART REFERENCE GUIDE IS AVAILABLE IN DVD FORMAT FOR YOUR CONVENIENCE.

THIS ITEM IS NON-RETURNABLE

THE COMPLETE ORTHODONTIC BONDING REFERENCE GUIDE

ORTHODONTIC BONDING REFERENCE GUIDE
DVD FORMAT

AD-1173

CHAIRSIDE BONDING TECHNIQUE CARDS



PRODUCT FEATURES:

INFORMATIVE AND CONVENIENT – ALL BONDING PROCEDURES AT YOUR FINGERTIPS! EVERY POSSIBLE INTRA-ORAL BONDING SITUATION IS COVERED. THESE DURABLE LAMINATED CARDS ARE ATTACHED ON A RING SO THEY CAN BE KEPT CHAIR SIDE.

CHAIRSIDE BONDING TECHNIQUE CARD SET

CHAIRSIDE BONDING TECHNIQUE CARDS

AD-1167

DISPOSABLE BENDA BRUSHES



1 1/2" BENDA BRUSHES (250)

AD-1172


SPONGE PELLETS



SPONGE PELLETS (500)

AD-1181

BRUSH HANDLE AND TIPS

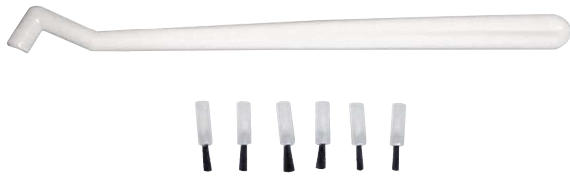

BRUSH HANDLE

AD-1170

BRUSH TIPS (250)

AD-1171

MIXING PADS



NON-ABSORBANT MIXING PADS (3)

AD-1175

NYLON MIXING SPATULAS



SPATULAS (50)

AD-1180

S.E.P. MIXING WELLS



MIXING WELLS WITH SHIELD (3)

AD-1023

MIXING WELLS



MIXING WELLS (5)

AD-1176

ALIGNER ATTACHMENT FORMS



OVAL SHAPE

AD-1026

RECTANGULAR SHAPE

AD-1027

ALIGNER ATTACHMENT KIT

AD-1040K

COMPOSITE AND CHAIN EXTRUSION KIT

AD-1168K